

Broad-range and Effective Detection of Human Noroviruses by Colloidal Gold Immunochromatographic Assay Based on the Shell Domain of the Major Capsid Protei

Meng Xu

Shanghai Jiao Tong University

Feifeng Lu

Shanghai Jiao Tong University

Qingping Wu

Guangdong Institute of Microbiology

Jumei Zhang

Guangdong Institute of Microbiology

Peng Tian

US Department of Agriculture

Liang Xue

Guangdong Institute of Microbiology

Chenang Lyu

Shanghai Jiao Tong University

Ting Xu

Shanghai Jiao Tong University

Dapeng Wang (✉ norovirus@163.com)

Shanghai Jiao Tong University

Research article

Keywords: human norovirus, broad-range detection, S domain protein, colloid gold, immunochromatographic assays

Posted Date: August 24th, 2020

DOI: <https://doi.org/10.21203/rs.3.rs-55525/v1>

License:   This work is licensed under a Creative Commons Attribution 4.0 International License.

[Read Full License](#)

Version of Record: A version of this preprint was published on January 11th, 2021. See the published version at <https://doi.org/10.1186/s12866-020-02084-z>.

Abstract

Background: Human noroviruses (HuNoVs) are a major cause of nonbacterial gastroenteritis in all ages worldwide. As the replication of HuNoVs *in vitro* is immature, the detection of HuNoVs is depended on molecular assays such as RT-PCR and RT-qPCR. However, these molecular-based techniques require special equipment, unique reagents, experienced operators to perform, and extended time to get results. In addition, the genotypes of HuNoVs are broad. Therefore, a method for rapidly, broad-range, and effective approach for HuNoVs would be more applicable for screening infected people when outbreaks occur.

Results: In this study, a colloidal-gold-based immunochromatographic assay (ICA) was developed for highly effective detection of HuNoVs in clinical samples. Monoclonal antibodies (McAbs) against the shell (S) domain in the major capsid protein of HuNoV were used in the ICA. The limitations of detection for HuNoVs in clinical samples were 1.2×10^6 genomic copies per gram (gc/g) and 4.4×10^5 gc/g for genogroup I and II (GI and GII) HuNoVs, respectively. A total of 122 samples were tested for HuNoVs by ICA and compared against that by RT-qPCR. The relative sensitivity, specificity and agreement of the ICA was 84.2 % (95% CI: 83.6-84.8 %), 100.0 % (95 % CI: 98.5-100.0 %) and 87.7 % (95% CI: 85.6-89.8 %), respectively. No cross-reaction with other common enteric viruses or bacteria was observed. The ICA could detect a broad range of genotypes, including GI.1, GI.3, GI.4, GI.6, GI.14, GII.2, GII.3, GII.4, GII.6, GII.13, and GII.17 HuNoVs.

Conclusions: Our results demonstrated that ICA targeting the S domain protein is a promising candidate for effectively improve identifying different genotypes of HuNoVs in clinical samples with high sensitivity and specificity.

Background

Human noroviruses (HuNoVs) are single-stranded RNA, non-enveloped viruses in the *Caliciviridae* family. The genome has three open reading frames (ORFs, 1 through 3), where ORF2 encodes a major capsid protein, VP1. VP1 consists of a shell (S) domain and a protruding (P) domain. The S protein, the most highly conserved region in VP1, forms a shell surrounding the RNA genome, while the P protein contains the most variable sequence [1-3]. HuNoVs are divided into 5 genogroups based on the VP1, including genogroup I (GI, n=9), GII (n=27), GIV (n=2), GVIII (n=1), and GIX (n=1) [4]. GI and GII HuNoVs are the epidemic strains circulating worldwide [5].

HuNoVs are the leading cause of epidemic gastroenteritis in the world [6-8]. The viral diseases are attributed to about 700 million infections and 200 thousand deaths globally, annually [9]. HuNoVs are highly-infectious, and infected patients can shed high titers of virus particles with strong resistance to the environmental factors [10, 11]. As cell culturing *in vitro* is still immature to replicate HuNoVs effectively [12], the detection of the virus mainly depends on the molecular approaches such as RT-PCR and RT-qPCR [13, 14]. These molecular approaches require specialized equipment, unique reagents, and significant time for sample preparation and assay execution, all of which make them impractical for either clinical

use or on-site assays of field samples requiring quick results. To efficiently control the spreading of HuNoVs in time, a highly effective, safe, and portable point-of-care assay is of great commercial and academic interest in monitoring the occurrence of outbreaks.

An immunochromatographic assay (ICA) is straightforward to screen samples for the presence of a target analyte, such as pathogens or biomarkers [15, 16]. In an ICA test, antibody against the viral capsid is labeled with gold particles. In the process of chromatography, the virus capsid is captured by another antibody coated on the solid-phase carrier, so that virus-antibody-colloidal gold particles combinations are also aggregated and show positive red results. ICAs are designed to be used without any equipment and require minimal training to do. The results can be read visually within a couple of minutes. Currently, several ICAs used for the detection of bacterial pathogens and viral pathogens have been reported [17-18]. However, many limitations of the ICA kit in HuNoVs detection were obvious, such as the detection of only specific genotypes of HuNoVs and the low sensitivity [19]. Therefore, a broad-range and highly effectively ICA was needed for point-of-care assay of HuNoVs.

Li et al. reported a McAb against S protein could cross-react with GI, GII, GIII, and GV of noroviruses (NoVs) [20]. Yoda et al. reported that two McAbs generated against the S domain of HuNoV (GII) capsid protein also recognized the viral protein of GI HuNoVs [21]. In addition, Parra et al. identified a broadly cross-reactive epitope in the S domain of the NoV capsid [2]. Therefore, the S domain of VP1, highly conserved, is a good candidate for detecting multiple strains of HuNoVs. However, few immunoassays detected the HuNoVs targeting the S protein. In this study, three McAbs were generated against the S domain of VP1 from a GII.4 HuNoV and used to develop an ICA for rapid and effective detection of both GI and GII HuNoVs.

Results

Selection of McAbs that recognize different epitopes on S proteins for ICA

The McAbs of S protein were named H9E, B4H, and J5D, further abbreviated as H, B, and J, respectively. The values of displacement factor (I) of B-H, H-B, B-J, J-B, H-J and J-H combinations were 32.8 %, 28.9 %, 34.1 %, 29.6 %, 31.2 % and 22.7 %, respectively. I value greater than 10.0 % indicates that two McAbs recognizes different epitopes on the S protein. Among the I values of six pairs, B-J had the best matching effect. The results were further tested by the ICA assay to find the best combination of McAbs binding to gold particles (Table 1). Finally, H was selected as the antibody for the labeled colloidal gold, and J as the T-line capture antibody.

Table 1 Results obtained from different combinations of the three antibodies on the colloidal gold platform

Colloidal-gold McAb	Test line McAb	Test line	Control line	
BH:	B4H	H9E	+	+
HB:	H9E	B4H	++	+
BJ:	B4H	J5D	+	++
JB:	J5D	B4H	++	+
HJ:	H9E	J5D	++	++
JH:	J5D	H9E	+	+

“+” represented the color depth. The score was evaluated by multiple operators with hidden labels.

Pretreatment methods to expose S proteins in the clinical samples

Two methods were tested for exposing the S domain in clinical samples. Under the temperature from 60 °C to 80 °C, the OD readings were gradually increased (roughly from 0.2 to 1.0). However, the OD readings were dramatically reduced when the temperature was greater than 90 °C (less than 0.2). A 3-minute at 80 °C treatment obtained the best result for heat-treatment (Fig. S1, see additional file 1). A progressive improvement of detection was displayed when pH was increased from 7.0 to 9.0 (optimal at pH=9.0, consistent with the best effect of heat treatment) (Fig. S2, see additional file 1). Better than heat-treatment, OD readings in alkali treatment samples were easy to control and stable over the time tested since there was no significant difference between the best 10-minute and the worst 5-minute treatment ($p = 0.131 > 0.05$). For simplicity and repeatability of the experimental operation, alkali treatment (samples were treated to pH=9.0 at room temperature for 10 min) was used to expose S proteins in all clinical samples tested by ICA.

The LOD of ICA

The C_t data and the copies were recorded in Table S1 (Additional file 2). The visual LOD of the ICA for the determination of S protein was 1.4 ng/ml. The LODs of viral genomic copies of sample 57404 (GI) and 1717 (GII) were 1.2×10^6 gc/g and 4.4×10^5 gc/g, respectively. The sensitivity test was performed in triplicate, and representative results were shown in Fig. 1.

Specificity of ICA

When clinical samples were detected by ICA, only samples containing HuNoVs showed positive results. Clinical samples containing other enteric viruses (four Rotaviruses, three Sapoviruses, two Astroviruses, and four Adenoviruses) or bacteria causing gastroenteritis (three *Salmonellas*) were all tested negative. Parts of the results were shown in Fig. 2.

Comparison of ICA with RT-qPCR

The comparison results of ICA and RT-qPCR were shown in Table 2. Eighty samples were positive by ICA from 95 RT-qPCR HuNoV positive samples. The sensitivity was 84.2 % (95% CI: 83.6-84.8 %) (80/95). No false-negative results indicated a specificity of 100.0 % (95% CI: 98.5-100.0 %) (27/27). The overall agreement of ICA and RT-qPCR was 87.7 % (95% CI: 85.6-89.8 %) (107/122). In addition, the titers of HuNoVs in 15 negative samples (1 GI.6, 1 GII.3, 12 GII.4, and 1 GII.6) were shown in Table S2 (Additional file 2). The viral loads of fifteen ICA false-negative samples were below the LOD of ICA.

Table 2 Comparison of HuNoVs detection in stool samples between the ICA kit and the RT-qPCR method

		ICA		
		Positive	Negative	Total
RT-qPCR	Positive	80	15	95
	Negative	0	27	27
Total		80	42	122

Sequencing results indicated that GII.4 was the dominant genotype (72/95) followed by GI.1 (4/95) in RT-qPCR HuNoV-positive samples. Other genotypes were GI.3 (2/95), GI.4 (3/95), GI.6 (2/95), GI.14 (1/95), GII.2 (2/95), GII.3 (3/95), GII.6 (3/95) GII.13 (1/95), and GII.17 (2/95) (Table 3). The ICA detected both GI and GII HuNoVs, including GI.1, GI.3, GI.4, GI.6, GI.14, GII.2, GII.3, GII.4, GII.6, GII.13, and GII.17 genotypes. The ICA detected both GI (11/12) and GII (69/83) HuNoVs including eleven genotypes.

Table 3 Detection results of genotypes by ICA

HuNoVs Genotypes	Positive Number in RT-qPCR	Number in ICA		Sensitivity (95 % CI)
		Positive	Negative	
GI.1	4	4	0	91.7 % (88.6-94.8)
GI.3	2	2	0	
GI.4	3	3	0	
GI.6	2	1	1	
GI.14	1	1	0	
(GI total)	(12)	(11)	(1)	
GII.2	2	2	0	83.1 % (80.0-86.3)
GII.3	3	2	1	
GII.4	72	60	12	
GII.6	3	2	1	
GII.13	1	1	0	
GII.17	2	2	0	
(GII total)	(83)	(69)	(14)	
Total	95	80	15	

Stability of ICA strips

Results demonstrated that acceptable activity (90.0 %) of the initial activity was retained for 21 days when the test strip was exposed at 60 °C. The most ability was equivalent to store at room temperature for approximately two years without apparent loss of activity.

Discussion

HuNoVs are one of the most influential foodborne pathogens in the world [22]. The burden of HuNoV outbreaks is high, reflected in the extensive infectivity and severe economic losses [23]. Acute gastroenteritis outbreaks associated with HuNoVs are particularly challenging to control because of their stability in the environment and efficient transmission in hospitals, hotels, schools, and domestic homes [24-26]. Therefore, a highly effective assay is needed to identify the HuNoVs for issuing emergency treatment.

Fecal specimens from patients could contain up to 10^9 gc/g of HuNoVs [27]. It is highly risky for the operators. Therefore, the inactivation of the viral particles in clinical samples is crucial before detection. Our results indicated that the capsid was greatly disassembled after treating by alkaline (pH 9.0 for 10

min) and heating (70.0 °C for 3 min). The alkali treatment was much better than heat treatment in terms of the stability of the results and the simplicity of operation. And the alkali environment also kept the colloidal gold in a stable state to bind McAbs.

The current molecular assays require special equipment, reagents, and experienced personnel. A simple approach with minimal hands-on time is more applicable and requisite as a screening assay. There were several commercial HuNoV detection kits available, such as Quick Ex-Norovirus, RIDA® QUICK (N1402), Immuno Search NV kit, NOROTOP+, and SD Bioline Norovirus [28-34]. However, these assays performed better for a single genotype. To detect multiple strains of HuNoVs, a couple of McAbs were selected for cross-reactions with various genotypes. A simple approach is urgently needed for broad-range genotypes of HuNoVs.

VP1 proteins (P plus S domains) could self-assemble to form virus-like particles (VLPs) which morphologically and antigenically resemble the viral particles [35, 36]. The protruding (P) domain of VP1 is highly stable and highly immunogenic similar to VLP [37]. The antigenicity of VLP and P particles is often strain-specific and is not suitable for virus screening. On the other hand, the S domain of the virus is highly conserved [7]. It was reported that two McAbs (1B4 and 1F6) of S protein against the capsid protein (NoV GI) were also found to react with NoV GI [21, 38]. In this study, McAbs against S protein were prepared to develop an ICA kit, which recognized VP1s from either GI or GII HuNoVs. The results demonstrated that our kit had good sensitivity and specificity, and detected both GI and GII HuNoVs (Table 3).

The epitopes on the S domain were hidden in the inner layer of viruses. They could not be recognized by the antibodies against S protein in the intact viral particles. It is necessary to expose S-protein for virus detection. In our ICA assay, the viral particles were disassembled by alkaline treatment and exposed S protein. It is an excellent way to minimize the risk of infection and much safer for testers than detection without inactivated samples. In addition, each viral particle contains 90 dimers of viral capsid [39]. After treatment, in theory, many more copies of S protein were exposed to the McAbs, and the sensitivity of the immunoassay could be improved. Therefore, it was a one-stone-two-bird solution.

It did report that IAC was used for the detection of NoV. IP-NoV kit and the kit developed by Takanashi was mainly for GII.3 and GII.4 HuNoVs [40-42]. Immuno Search NV kit was developed for several genotypes (GI.1, GI.11, GII.2, GII.3, GII.4, GII.5, and GII.6 VLPs), but the limit of detection was not shown [29]. SD Bioline Norovirus could detect GI.3 for GI but except GII.2, GII.6, and GII.16 for GII HuNoVs [31, 32]. Limited genotypes of HuNoVs were detected by those kits. It is noteworthy that our ICA could detect a much broader range of genotypes (including GI.1, GI.3, GI.4, GI.6, GI.14, GII.2, GII.3, GII.4, GII.6, GII.13 and GII.17) (Table 3).

Compared to RT-qPCR, the sensitivity, specificity and agreement of the ICA were 84.2 % (95% CI: 83.6-84.8 %), 100.0 % (95% CI: 98.5-100.0 %) and 87.7 % (95% CI: 85.6-89.8 %), respectively (Table 2). The overall results were competitive with other commercial kits. The sensitivity (84.2%) of our ICAs was higher than IP-NoV kit (ranged from 72.7 % to 78.9 %), Quick Ex-Norovirus (54.5 %), RIDA® QUICK (from 68.8 % to

82.5 %), Immuno Search NV kit (75.4 %) and NOROTOP+ (51.4 %) (See Table S3, Additional file 3). In addition, Rotavirus, Sapovirus, Astrovirus, Adenovirus, *Salmonella*, and *Vibrio parahaemolyticus* are the most common non-norovirus viral and bacterial pathogens that cause gastroenteritis in China [43]. None of the 16 specimens were tested positive in the evaluation of specificity in our study.

These performance characteristics of ICA kits were dependent on the tested individual and the titer of HuNoVs in the sample. There were several possibilities for false-negative results, including lower viral titer in the tested clinical samples and sensitivity of the McAbs selected in rapid ICA tests. In our study, 15 false-negative samples had viral titers below the LOD. In the future, the LOD of our kit should be improved by quantum dot or signal amplification system. We believe the ICA is a promising approach for screening HuNoVs from clinical samples, especially treating the outbreaks.

Conclusions

In this study, an ICA targeting the S-domain of the HuNoV VP1 protein for the detection of HuNoVs was developed. This assay was able to reduce detection times and achieve a sensitivity of 84.2 % (95% CI: 83.6-84.8 %) and a specificity of 100.0 % (95% CI: 98.5-100.0 %), which is competitive to other immunological kits without using special equipment and reagents. We believe this assay represents a promising approach for screening HuNoVs from the field, clinical samples or both.

Methods

Preparation of VP1s and S proteins

HuNoV GI.1 and GII.4 ORF2 gene fragments (GenBank Nos. M87661 and KM114291, respectively) were subcloned separately into the pET-28a prokaryotic expression plasmid. HuNoV GI and GII VP1 proteins were expressed from respective recombinant plasmid vectors pET28a-ORF2 GI.1 and GII.4 in *Escherichia coli* (*E. coli*) BL21 cells [44].

HuNoV GII.4 ORF2 (GenBank No. KM114291) was used as a template to amplify the nucleic acid fragment of the S domain. The upstream primer was 5'-GAATTCATGAAGATGGCGTCGAGTG-3'; the downstream primer was 5'-CTCGAGCTCAACTGTGGGTGGCAC-3'. *EcoR* I and *Xho* I restriction sites were appended to the 5' ends of forward and reverse primers, respectively. After amplification and digestion, the nucleic acid fragment of the S domain was inserted into a pSmart vector (Frdbio, China) to generate recombinant pSmart-S. The recombinant S protein was induced and expressed, as described in the previous report [45].

Preparation of anti-S protein McAbs

The S domain McAbs were prepared as we previously described [46]. Three McAbs (H9E, B4H, and J5D, further abbreviated as H, B, and J, respectively) were selected for their ability to bind pSmart-S expression

product, but not with *E. coli* residual proteins. These three McAbs also demonstrated specificity against both GI and GII HuNoV VP1s via Western Blot.

Selection of McAbs that recognize different epitopes of S protein by ELISA assay

Each reaction of the ICA assay utilizes a labeled McAb and another immobilized (“capture”) McAb. All six pair permutations of the three McAbs (B-H, H-B, B-J, J-B, H-J, and J-H) were tested. First, the three McAbs (B, H, and J) were two-fold serially diluted and tested for their single epitope saturation concentrations using recombinant S protein (100.0 µg/well). A concentration one step lower than the concentration with a significant decrease in OD₄₅₀ value was defined as a single epitope saturation concentration and named as B₁ (or H₁, J₁). Then S protein (100.0 µg/well) was coated overnight at 4 °C, and a 100.0 µl epitope-saturated McAb (e.g., B) was incubated at 37 °C for 1 h. After washing, 100.0 µl second saturated McAb (e.g., H) in the pair tested was incubated under the same conditions. OD₄₅₀ was recorded as BH₂. Displacement factor (*I*) was defined as the ratio of the increased effect of cross-reactivity over the separate effect of the second added antibody. If $I > 10\%$ [47], it indicates that the recognition sites of the two McAbs are distinct. *I* values were calculated according to the formula $I_{(e.g., BH)} = (BH_2 - B_1)/H_1 \times 100\%$.

Selection of McAbs used on the colloidal gold platform

The performance of ICA depended critically upon the combination of the optimal antibody sandwich pair with colloidal gold. ELISA results needed to be considered in conjunction with ICA testing. The capture antibody on the Test line (T line) and the control antibody (goat-anti-mouse Ig G, Beyotime Biotech, Shanghai, China) on the Control line (C line) (Fig. 3) was both 2.0 mg/ml, and the labeled antibody was exceeded (40 µg/ml). The concentration of S protein used in the test was 0.3 µg/ml. The final selection was based on the combination of the color effect of both C and T lines.

Source of HuNoVs clinical samples (including suspected HuNoVs)

A total of 122 fecal specimens were collected and tested. Five HuNoVs clinical samples [57404 (GI.1), 3010 (GI.1), 1704 (GII.4), 1717 (GII.4), 1028 (GII.4)] were kindly provided by Dr. Ningbo Liao (Zhejiang Provincial Center for Disease Control and Prevention (CDC)). Total 117 clinical diarrheal samples which kindly provided by the Affiliated Hospital of Guangzhou Medical University (N=3), the Affiliated Hospital of Sun Yat-sen University (N=7), Chinese CDC (N=39) and Anhui Provincial CDC (N=68) were also tested by both RT-qPCR and ICA. Clinical samples provided by Zhejiang Provincial and Chinese CDC were from patients with acute gastroenteritis in 2018. The samples of Anhui provincial CDC and Affiliated Hospital of Guangzhou Medical University were all in 2017. The samples of the Affiliated Hospital of Sun Yat-sen University were from 2015 to 2017.

Detection of HuNoVs by RT-qPCR and calculation of viral genomic copies

Real-time RT-PCR was performed using a commercial one-step RT-qPCR kit (Sangon Biotech, China), consisting of 12.5 µl 2 × one-step RT-qPCR Master Mix (with SYBR Green), 0.65 µl RT enzyme Mix, 0.4 µl

for each primer (0.16 $\mu\text{mol/l}$) (See Table S4, Additional file 4) [48]. 2.0 μl RNA template and RNase free ddH₂O to a total volume of 25.0 μl . Amplification reactions were as follows: reverse transcription at 50 °C for 30 min; heat-denatured at 95 °C for 3 min; 40 cycles with denaturation at 95 °C for 10 s, annealing, and extension at 60 °C for 30 s. Fluorescence signals were collected at the end of each extension step. The highest dilution of real-time RT-PCR to generate a positive cycle threshold (C_t) signal was one real-time RT-PCR unit (RT-qPCR) [13]. C_t values versus \log_{10} viral genomic copies linear standard curves were generated from a continuous 10-fold dilution (See Figure S3, Additional file 5). Amplified DNA fragments were sequenced on ABI 3730XL at Personalbio (Shanghai, China). Automated genotypes were analyzed with Norovirus Typing Tool Version 2.0 (www.rivm.nl/mpf/norovirus/typingtool).

Conditions for exposing S protein from viral capsid in clinical samples

Pretreatment comprises of physical or chemical methods. Heat-denaturation was used for the physical treatment. Briefly, ten percent (w/v) fecal sample diluted with ddH₂O was put in a water bath at a set temperature (60 °C, 70 °C, 80 °C, 90 °C, 100 °C) for 1 to 5 min with an interval of 1 min. The reducing agent (DTT) (Thermo Fisher, China) and alkaline were used for chemical treatment. The pH in PBS was adjusted to 6.0 to 10.0. Fecal samples were dissolved to 10.0 % with PBS buffers (pH 6.0, 7.0, 8.0, 9.0 or 10.0) at room temperature (25 °C) for 5, 10, 15, 20, 25 and 30 minutes. DTT was added to the homogenized solution of clinical sample at a final concentration of 1.0 %. A sandwich ELISA was undertaken to detect HuNoVs, as described previously [49].

Performing an ICA test

The procedure of preparing ICA strips was described in the Supplementary Material (Figure S4, Table S5-7, See Additional file 6). The sample preparation process was as follows: 10 % (w/v) homogenized solution of stool slurries was prepared with PBS (pH 9.0) (10,000 \times g, 10 min). DTT was added to a final concentration (1.0 %). After 10 minutes at 25 °C, 50 μl mixtures were added to the sample pad. If enough HuNoV antigens were present, the binding of the gold-labeled antigen complex occurred at both the T and C lines. The presence of C line confirmed that the test was valid.

Limit of detection (LOD) of ICA for S protein and clinical samples

Purified S protein was used at a concentration of 22.4 ng/ml, 11.2 ng/ml, 5.6 ng/ml, 2.8 ng/ml, 1.4 ng/ml and 0.7 ng/ml. PBS buffer was used as blank control. Two HuNoV positive samples (57404 GI.1 and 1717 GII.4) were two-fold diluted with a range of 1.6×10^5 to 5×10^6 gc/g (GI) and 1.1×10^5 to 3.5×10^6 gc/g (GII). All analyses were conducted in triplicate to determine the repeatability of the results.

Evaluation of the specificity of ICA

To further evaluate the specificity of antibodies in the colloidal gold test, four Rotavirus, three Sapovirus, two Astrovirus and four Adenovirus (stool sample provided by Zhejiang CDC, the affiliated Hospital of Guangzhou Medical University and the affiliated Hospital of Sun Yat-sen University) and three *Salmonella*

(ATCC 14028, CMCC 50115 and CICC 21482) were tested. The copies of Rotavirus, Sapovirus, Astrovirus, and Adenovirus RNA were more than 10^7 gc/g feces.

Stability test

To determine the stability of ICA strips, they were examined by thermal acceleration tests at 60 °C. The activity of the antibodies on each assay was tested by using the lowest detectable S protein concentration (1.4 ng/ml). The identical strips were tested at appropriate intervals (1, 2, 3, 7, 14, 21, and 28 days). By comparing the result of Sylvia's report [50], the shelf-life could be estimated at room temperature.

Statistical analysis

One-way ANOVA was used for data analysis. SPSSAU, an online data analysis tool, was used for statistical analysis (www.spssau.com).

Abbreviations

CI: confidence interval; ATCC: American type culture collection; ELISA: enzyme-linked immunosorbent assay; DTT: dithiothreitol

Declarations

Acknowledgments

We would extend special thanks to Dr. Ningbo Liao from Zhejiang Provincial CDC, Dr. Miao Jin from Chinese CDC, and Dr. Yong Sun from Anhui Provincial CDC for kindly providing test samples. We also greatly appreciated David Yang (PSMRU, WRRRC, ARS, USDA) for his assistance in revising the manuscript and Weicheng Cai (Guangdong Institute of Microbiology) for his complete support of the experiment.

Funding

This work was jointly supported by the Ministry of Science and Technology of the People's Republic of China (2017YFC1601200), the National Natural Science Foundation of China (31772078) and the Agri-X Interdisciplinary Fund of Shanghai Jiao Tong University (Agri-X 2017005).

Availability of data and materials

The datasets supporting the conclusions of this article are included within the article and its additional files. The datasets used and/or analyzed during the current study are available from the corresponding author on reasonable request.

Authors' contributions

MX: methodology, formal analysis, and writing the original draft. FL: software, investigation, and validation. QW: methodology and supervision. JZ: project administration and data curation. PT: conceptualization, reviewing, editing, and validation. LX: resources and supervision. CL: resources and visualization. TX: resources and visualization. DW: conceptualization, supervision, and funding acquisition. All authors read and approved the final manuscript.

Ethics approval and consent to participate

Not applicable

Consent for publication

Not applicable

Competing interests

The authors declare that they have no conflict of interest.

References

1. Parker TD, Kitamoto N, Tanaka T, Hutson AM, Estes MK. Identification of genogroup I and genogroup II broadly reactive epitopes on the norovirus capsid. *J Virol*. 2005;79(12):7402-9. doi: 10.1128/JVI.79.12.7402-7409.2005
2. Parra GI, Azure J, Fischer R, Bok K, Sandoval-Jaime C, Sosnovtsev SV, Sander P, Green KY. Identification of a broadly cross-reactive epitope in the inner shell of the Norovirus capsid. *PLoS One*. 2013;8(6):e67592. doi: 10.1371/journal.pone.0067592
3. Parra GI, Sosnovtsev SV, Abente EJ, Sandoval-Jaime C, Bok K, Dolan MA, Green KY. Mapping and modeling of a strain-specific epitope in the Norwalk virus capsid inner shell. *Virology*. 2016;492:232-41. doi: 10.1016/j.virol.2016.02.019
4. Chhabra P, de Graaf M, Parra GI, Chan MC, Green K, Martella V, Wang Q, White PA, Katayama K, Vennema H, Koopmans MPG, Vinjé J. Updated classification of norovirus genogroups and genotypes. *J Gen Virol*. 2019;100(10):1393-1406. doi: 10.1099/jgv.0.001318
5. Motoya T, Nagasawa K, Matsushima Y, Nagata N, Ryo A, Sekizuka T, Yamashita A, Kuroda M, Morita Y, Suzuki Y, Sasaki N, Katayama K, Kimura H. Molecular evolution of the VP1 gene in human norovirus GII.4 variants in 1974-2015. *Front Microbiol*. 2017;8:2399. doi: 10.3389/fmicb.2017.02399
6. Glass RI, Parashar UD, Estes MK. Norovirus gastroenteritis. *N Engl J Med*. 2009;361(18):1776–1785. doi:10.1056/NEJMr0804575
7. Kabue JP, Meader E, Hunter PR, Potgieter N. Human norovirus prevalence in Africa: a review of studies from 1990 to 2013. *Trop Med Int Health*. 2016;21(1):2-17. doi: 10.1111/tmi.12617
8. Kowalzik F, Binder H, Zöller D, Riera-Montes M, Clemens R, Verstraeten T, Zepp F. Norovirus gastroenteritis among hospitalized patients, Germany, 2007-2012. *Emerg Infect Dis*.

- 2018;24(11):2021-2028. doi: 10.3201/eid2411.170820
9. Kirk MD, Pires SM, Black RE, Caipo M, Crump JA, Devleeschauwer B, Döpfer D, Fazil A, Fischer-Walker CL, Hald T, Hall AJ, Keddy KH, Lake RJ, Lanata CF, Torgerson PR, Havelaar AH, Angulo FJ. World health organization estimates of the global and regional disease burden of 22 foodborne bacterial, protozoal, and viral diseases, 2010: A Data Synthesis. *PLoS Med.* 2015;12(12):e1001921. doi: 10.1371/journal.pmed.1001921
 10. de Graaf M, van Beek J, Koopmans MP. Human norovirus transmission and evolution in a changing world. *Nat Rev Microbiol.* 2016;14(7):421-33. doi:10.1038/nrmicro.2016.48
 11. Teunis PF, Moe CL, Liu P, Miller SE, Lindesmith L, Baric RS, Le Pendu J, Calderon RL. Norwalk virus: how infectious is it? *J Med Virol.* 2008;80(8):1468-76. doi: 10.1002/jmv.21237
 12. Costantini V, Morantz EK, Browne H, Ettayebi K, Zeng XL, Atmar RL, Estes MK, Vinjé J. Human norovirus replication in human intestinal enteroids as model to evaluate virus inactivation. *Emerg Infect Dis.* 2018;24(8):1453-1464. doi: 10.3201/eid2408.180126
 13. Radin D, D'Souza DH. Simple and rapid detection of Human norovirus from produce using SYBR Green I-based real-time RT-PCR. *Food Environ Virol* 2011;3,121–129. doi: 10.1007/s12560-011-9066-5
 14. Zhou Z, Tian Z, Li Q, Tian P, Wu Q, Wang D, Shi X. *In situ* capture RT-qPCR: A new simple and sensitive method to detect human norovirus in oysters. *Front Microbiol.* 2017;8:554. doi: 10.3389/fmicb.2017.00554
 15. Nakayama T, Zhao J, Takeuchi D, Kerdsin A, Chiranairadul P, Areeratana P, Loetthong P, Pienpringam A, Akeda Y, Oishi K. Colloidal gold-based immunochromatographic strip test compromising optimised combinations of anti-*S. suis* capsular polysaccharide polyclonal antibodies for detection of *Streptococcus suis*. *Biosens Bioelectron.* 2014;60:175-9. doi: 10.1016/j.bios.2014.03.074
 16. Yu G, Yu X, Yang G, Tang Y, Diao Y. A novel diagnostic method to detect Duck Tembusu virus: a colloidal gold-based immunochromatographic assay. *Front Microbiol.* 2018;9:1001. doi: 10.3389/fmicb.2018.01001
 17. Huang SH, Wei HC, Lee YC. One-step immunochromatographic assay for the detection of *Staphylococcus aureus*. *Food Control.* 2007;18(8):893-897. doi: 10.1016/j.foodcont.2006.05.005
 18. Souza IO, Schaer RE, de M. Lima FW, Raynal JT, Sousa TJ, Silva AA, Sardi SI, Meyer R, Portela RW. Validation of an immunochromatographic assay for the multiple detection of specific antibodies against HIV, HBV, and HCV. *Clin Vaccine Immunol.* 2013;20(1):128. doi: 10.1128/CVI.00485-12
 19. Kirby A, Gurgel RQ, Dove W, Vieira SC, Cunliffe NA, Cuevas LE. An evaluation of the RIDASCREEN and IDEIA enzyme immunoassays and the RIDAQUICK immunochromatographic test for the detection of norovirus in faecal specimens. *J Clin Virol.* 2010;49(4):254-7. doi: 10.1016/j.jcv.2010.08.004.
 20. Li X, Zhou R, Tian X, Li H, Zhou Z. Characterization of a cross-reactive monoclonal antibody against Norovirus genogroups I, II, III and V. *Virus Res.* 2010;151(2):142-7. doi: 10.1016/j.virusres.2010.04.005

21. Yoda T, Suzuki Y, Terano Y, Yamazaki K, Sakon N, Kuzuguchi T, Oda H, Tsukamoto T. Precise characterization of norovirus (Norwalk-like virus)-specific monoclonal antibodies with broad reactivity. *J Clin Microbiol.* 2003;41(6):2367-71. doi: 10.1128/jcm.41.6.2367-2371.2003
22. Dewey-Mattia D, Manikonda K, Hall AJ, Wise ME, Crowe SJ. Surveillance for foodborne disease outbreaks - United States, 2009-2015. *MMWR Surveill Summ.* 2018;67(10):1-11. doi: 10.15585/mmwr.ss6710a1
23. Verstraeten T, Cattaert T, Harris J, Lopman B, Tam CC, Ferreira G. Estimating the burden of medically attended Norovirus gastroenteritis: modeling linked primary care and hospitalization datasets. *J Infect Dis.* 2017;216(8):957-965. doi: 10.1093/infdis/jix410
24. Wikswo ME, Hall AJ; Outbreaks of acute gastroenteritis transmitted by person-to-person contact—United States, 2009-2010. *MMWR Surveill Summ.* 2012;61(9):1-12. doi: [10.2105/AJPH.2014.10411e13](https://doi.org/10.2105/AJPH.2014.10411e13)
25. Atmar RL, Opekun AR, Gilger MA, Estes MK, Crawford SE, Neill FH, Graham DY. Norwalk virus shedding after experimental human infection. *Emerg Infect Dis.* 2008;14(10):1553-7. doi: 10.3201/eid1410.080117
26. Johnston CP, Qiu H, Ticehurst JR, Dickson C, Rosenbaum P, Lawson P, Stokes AB, Lowenstein CJ, Kaminsky M, Cosgrove SE, Green KY, Perl TM. Outbreak management and implications of a nosocomial norovirus outbreak. *Clin Infect Dis.* 2007;45(5):534-40. doi: [10.1086/520666](https://doi.org/10.1086/520666)
27. Buckow R, Isbarn S, Knorr D, Heinz V, Lehmacher A. Predictive model for inactivation of feline calicivirus, a norovirus surrogate, by heat and high hydrostatic pressure. *Appl Environ Microbiol.* 2008;74(4):1030-8. doi: 10.1128/AEM.01784-07
28. Mutoh K, Hakamata A, Yagi H, Kurokawa K, Miki N, Kurita I. Evaluation of new commercial immunochromatography kit for norovirus in feces. *Pediatr Int.* 2009;51(1):164-6. doi: 10.1111/j.1442-200X.2008.02788.x.
29. Bruggink LD, Dunbar NL, Marshall JA. Evaluation of the updated RIDA®QUICK (Version N1402) immunochromatographic assay for the detection of norovirus in clinical specimens. *J Virol Methods.* 2015;223:82-7. doi: 10.1016/j.jviromet.2015.07.019.
30. Khamrin P, Takanashi S, Chan-It W, Kobayashi M, Nishimura S, Katsumata N, Okitsu S, Maneekarn N, Nishio O, Ushijima H. Immunochromatography test for rapid detection of norovirus in fecal specimens. *J Virol Methods.* 2009;157(2):219-22. doi: 10.1016/j.jviromet.2008.12.012
31. Ambert-Balay K, Pothier P. Evaluation of 4 immunochromatographic tests for rapid detection of norovirus in faecal samples. *J Clin Virol.* 2013;56(3):194-8. doi: 10.1016/j.jcv.2012.11.001
32. Kim HS, Hyun J, Kim JS, Song W, Kang HJ, Lee KM. Evaluation of the SD Biotec Norovirus rapid immunochromatography test using fecal specimens from Korean gastroenteritis patients. *J Virol Methods.* 2012;186(1-2):94-8. doi: 10.1016/j.jviromet.2012.08.014
33. Park KS, Baek KA, Kim DU, Kwon KS, Bing SH, Park JS, Nam HS, Lee SH, Choi YJ. Evaluation of a new immunochromatographic assay kit for the rapid detection of norovirus in fecal specimens. *Ann Lab Med.* 2012;32(1):79-81. doi: 10.3343/alm.2012.32.1.79

34. Geginat G, Kaiser D, Schrempf S. Evaluation of third-generation ELISA and a rapid immunochromatographic assay for the detection of norovirus infection in fecal samples from inpatients of a German tertiary care hospital. *Eur J Clin Microbiol Infect Dis*. 2012;31(5):733-7. doi: 10.1007/s10096-011-1366-z
35. Tan M, Hegde RS, Jiang X. The P domain of norovirus capsid protein forms dimer and binds to histo-blood group antigen receptors. *J Virol*. 2004;78(12):6233-42. doi: 10.1128/JVI.78.12.6233-6242.2004
36. Jiang X, Wang M, Graham DY, Estes MK. Expression, self-assembly, and antigenicity of the Norwalk virus capsid protein. *J Virol*. 1992;66(11):6527-32. doi: 10.0000/PMID1328679
37. Tan M, Jiang X. Norovirus P particle: a subviral nanoparticle for vaccine development against norovirus, rotavirus and influenza virus. *Nanomedicine (Lond)*. 2012;7(6):889-97. doi: 10.2217/nnm.12.62
38. Yoda T, Terano Y, Shimada A, Suzuki Y, Yamazaki K, Sakon N, Oishi I, Utagawa ET, Okuno Y, Shibata T. Expression of recombinant Norwalk-like virus capsid proteins using a bacterial system and the development of its immunologic detection. *J Med Virol*. 2000;60(4):475-81. doi: [1002/\(SICI\)1096-9071\(200004\)60:4<475::AID-JMV17>3.0.CO;2-B](https://doi.org/10.1002/(SICI)1096-9071(200004)60:4<475::AID-JMV17>3.0.CO;2-B)
39. Prasad BV, Hardy ME, Dokland T, Bella J, Rossmann MG, Estes MK. X-ray crystallographic structure of the Norwalk virus capsid. *Science*. 1999;286(5438):287-90. doi: 10.1126/science.286.5438.287
40. Nguyen TA, Khamrin P, Takanashi S, Le Hoang P, Pham le D, Hoang KT, Satou K, Masuoka Y, Okitsu S, Ushijima H. Evaluation of immunochromatography tests for detection of rotavirus and norovirus among Vietnamese children with acute gastroenteritis and the emergence of a novel norovirus GII.4 variant. *J Trop Pediatr*. 2007;53(4):264-9. doi: 10.1093/tropej/fmm021
41. Thongprachum A, Khamrin P, Chaimongkol N, Malasao R, Okitsu S, Mizuguchi M, Maneekarn N, Ushijima H. Evaluation of an immunochromatography method for rapid detection of noroviruses in clinical specimens in Thailand. *J Med Virol*. 2010;82(12):2106-9. doi: 10.1002/jmv.21916
42. Takanashi S, Okame M, Shiota T, Takagi M, Yagyu F, Tung PG, Nishimura S, Katsumata N, Igarashi T, Okitsu S, Ushijima H. Development of a rapid immunochromatographic test for noroviruses genogroups I and II. *J Virol Methods*. 2008;148(1-2):1-8. doi: 10.1016/j.jviromet.2007.10.010
43. Xue Y, Pan H, Hu J, Wu H, Li J, Xiao W, Zhang X, Yuan Z, Wu F. Epidemiology of norovirus infections among diarrhea outpatients in a diarrhea surveillance system in Shanghai, China: a cross-sectional study. *BMC Infect Dis*. 2015;15:183. doi: 10.1186/s12879-015-0922-z
44. Niu M, Yu Q, Tian P, Gao Z, Wang D, Shi X. Engineering bacterial surface displayed Human norovirus capsid proteins: a novel system to explore interaction between Norovirus and ligands. *Front Microbiol*. 2015;6:1448. doi: 10.3389/fmicb.2015.01448
45. Xu Q, Ni P, Liu D, Yin Y, Li Q, Zhang J, Wu Q, Tian P, Shi X, Wang D. A bacterial surface display system expressing cleavable capsid proteins of Human norovirus: a novel system to discover candidate receptors. *Front Microbiol*. 2017;8:2405. doi: 10.3389/fmicb.2017.02405
46. Wang D, Wu Q, Kou X, Yao L, Zhang J. Distribution of norovirus in oyster tissues. *J Appl Microbiol*. 2008;105(6):1966-72. doi: 10.1111/j.1365-2672.2008.03970.x

47. Xu Y, van Huystee R. Identification of an antigenic determinant on anionic peanut peroxidase by monoclonal antibodies. *J Exp Bot.* 1991;42(7):935-945. doi: 10.1093/jxb/42.7.935
48. Kageyama T, Kojima S, Shinohara M, Uchida K, Fukushi S, Hoshino FB, Takeda N, Katayama K. Broadly reactive and highly sensitive assay for Norwalk-like viruses based on real-time quantitative reverse transcription-PCR. *J Clin Microbiol.* 2003;41(4):1548-57. doi: 10.1128/jcm.41.4.1548-1557.2003
49. Yoda T, Terano Y, Suzuki Y, Yamazaki K, Oishi I, Kuzuguchi T, Kawamoto H, Utagawa E, Takino K, Oda H, Shibata T. Characterization of Norwalk virus GI specific monoclonal antibodies generated against *Escherichia coli* expressed capsid protein and the reactivity of two broadly reactive monoclonal antibodies generated against GII capsid towards GI recombinant fragments. *BMC Microbiol.* 2001;1:24. doi: 10.1186/1471-2180-1-24
50. Sylvia PO, Carmen LS, Montserrat RU. Shelf-life prediction of an infant formula using an accelerated stability test (Rancimat). *J. Agr. Food Chem.* 1995;43(11):2879-2882. doi: 10.1021/jf00059a020

Figures



Figure 1

LOD of ICA for S protein and clinical samples (A) Sensitivity for S proteins (twofold dilutions from 22.4 ng/ml to 0.7 ng/ml), (B) Sample 57404 (GI.1) (twofold dilutions from 5.0×10^6 to 1.6×10^5 gc/g) and (C) Sample 1717 (GII.4) of different virus copies (twofold dilutions from 3.5×10^6 to 1.1×10^5 gc/g) were detected with test strips. PBS elution buffer (pH 7.4) was used as blank control.

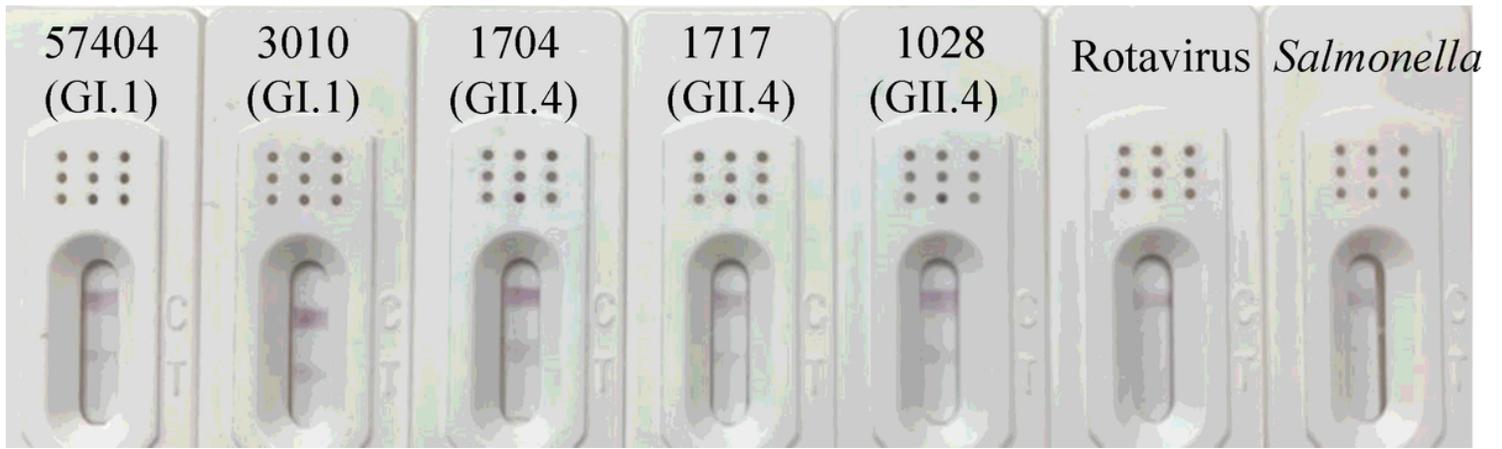


Figure 2

Specificity of ICA Five HuNoVs clinical samples, Rotavirus, and Salmonella cultured samples were tested with the strips.

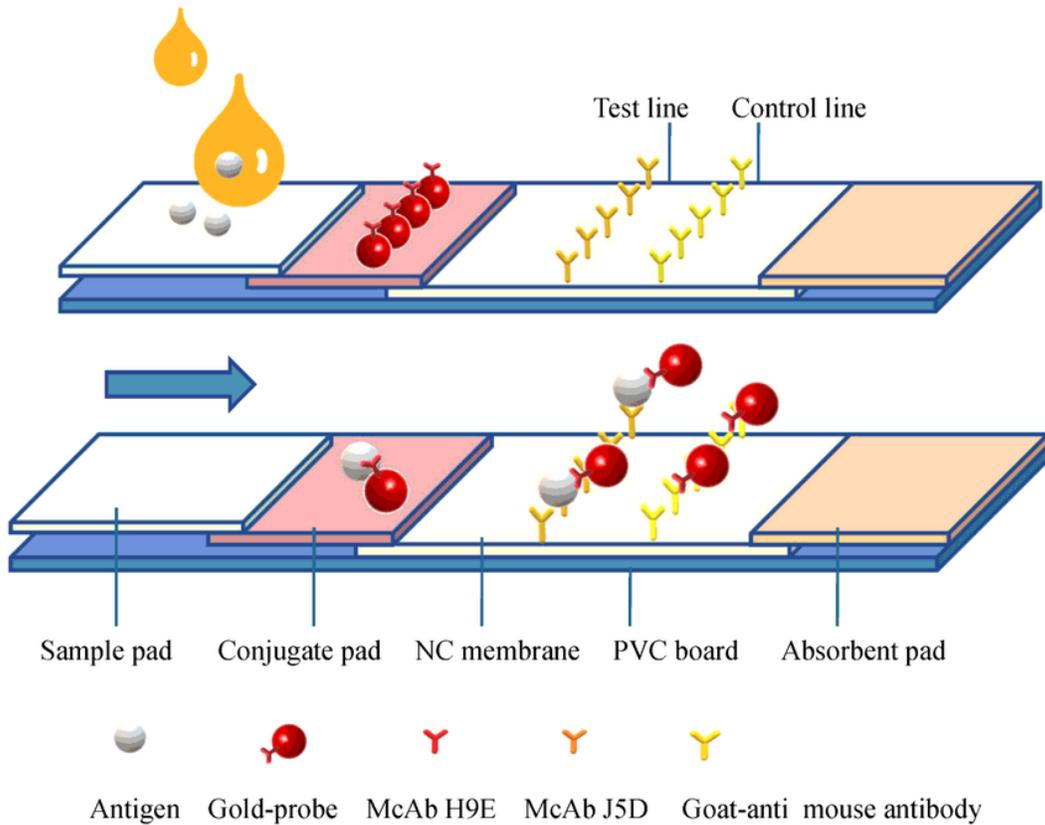


Figure 3

Schematic diagram of ICA test. McAb H9E was used as labeled antibodies with colloidal gold particles (in red); McAb J5D was used as capture-antibodies in the T line; goat-anti mouse McAb Ig G was used in the C line. Arrow indicates the direction of the movement of antigens (capsid protein of noroviruses).

Supplementary Files

This is a list of supplementary files associated with this preprint. Click to download.

- [Additionalfile6Fig.S4TableS57.docx](#)
- [Additionalfile5Fig.S3.docx](#)
- [Additionalfile4TableS4.docx](#)
- [Additionalfile3TableS3.docx](#)
- [Additionalfile2TableS12.docx](#)
- [Additionalfile1Fig.S12.docx](#)